

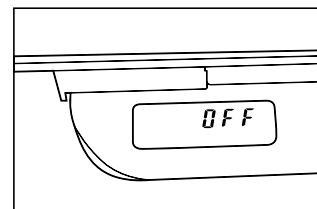
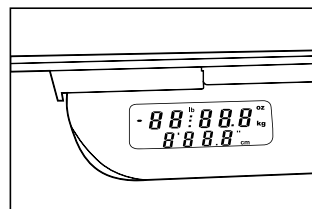
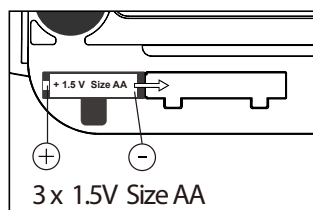



## Installation

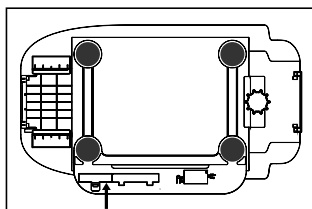

## Turn on

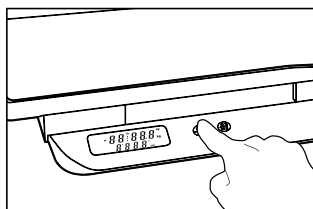


OR

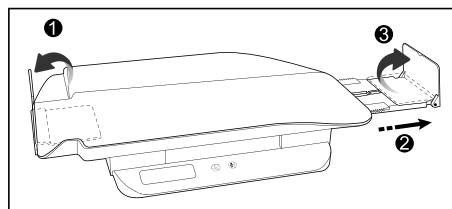


Auto-turn on only works when ruler is fully extended

## Turn off

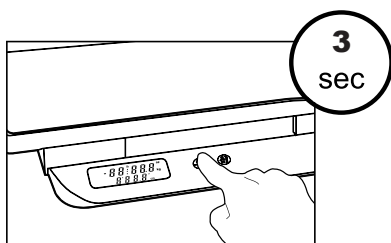


OR

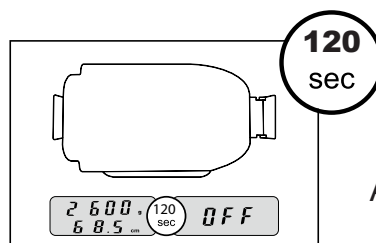


Auto-turn off after 2 minutes


## Measure weight + height

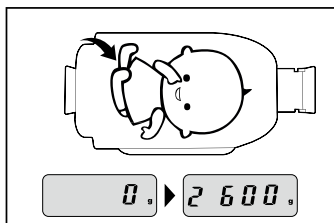

Put baby on scale to measure weight

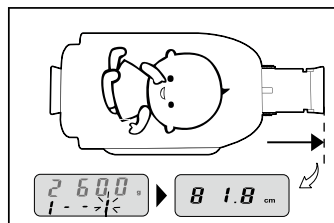


Fully extend the digital ruler

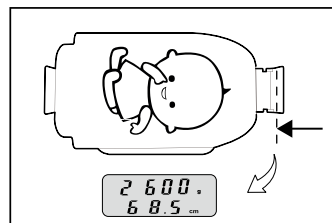


Extend the flap and push in the ruler until it touches your baby's head


## Measure height only

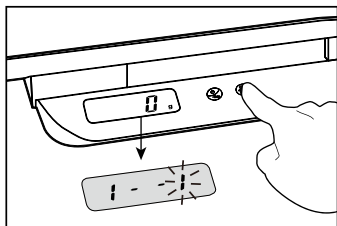

Press once

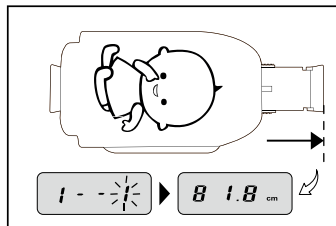


Fully extend the digital ruler


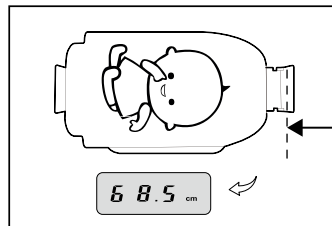


Extend the flap and push in the ruler until it touches your baby's head


## Tare (set to 0 even with cloth on the scale)

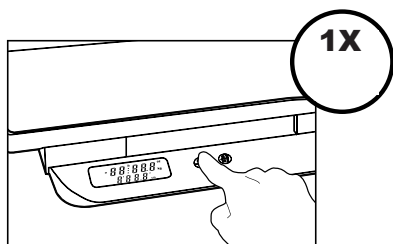


## Switch between grams and ounce

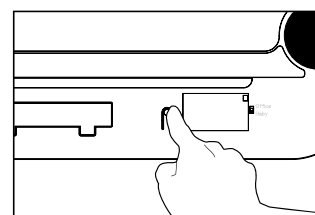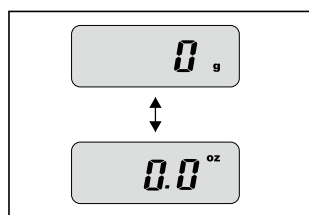

## Messages

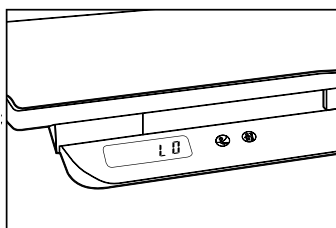

Low battery  
Replace the batteries

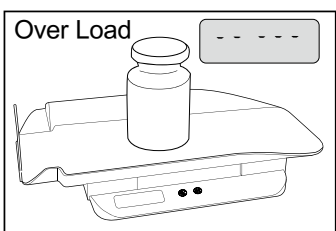

Overload  
Maximum 20 kg can be measured

## SPECIFICATIONS

Power supply:..... 3x 1,5Volt, size AA  
Battery lifetime:..... > 1 year  
Maximum load:..... 20 Kg  
Gradation:..... 5 grams  
Maximum length: ..... 81,8 cm  
Length gradation: ..... 0,2 cm  
Work temperature:..... 10 - 40°C

Service

Help

↓

[WWW.ALECTO.NL](http://WWW.ALECTO.NL)  
[SERVICE@ALECTO.NL](mailto:SERVICE@ALECTO.NL)

Hesdo, Australiëlaan 1  
5232 BB, 's-Hertogenbosch  
The Netherlands